


ACUPRESSURE LUMBAR BELT WITH VIBRATION AND MASSAGE

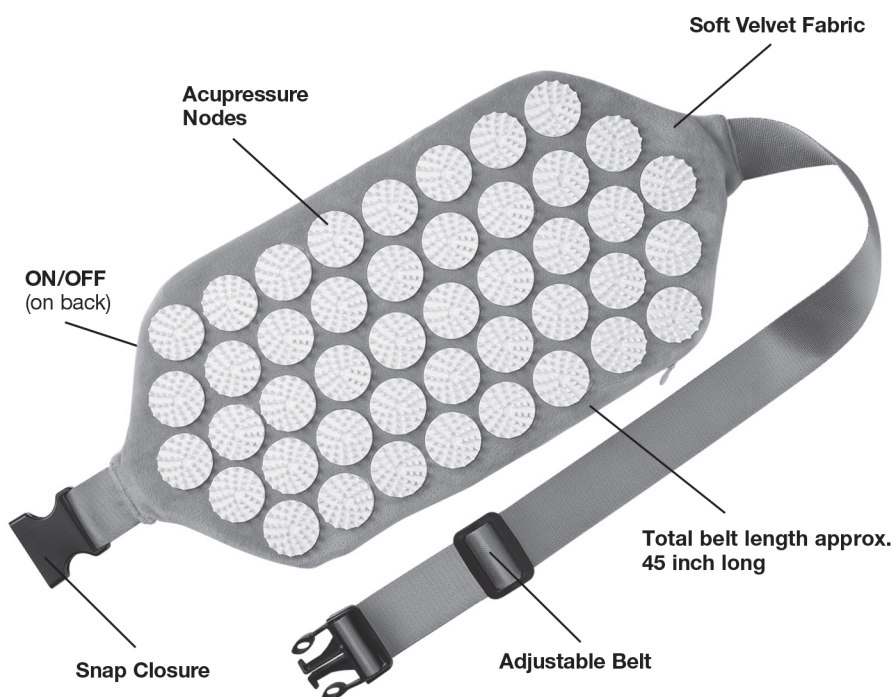

Powerful and Portable
for home, office or car

Battery Operated
requires 2 AA batteries,
not included

ONE YEAR LIMITED WARRANTY

Conair will (at our option) repair or replace your appliance free of charge for 12 months from the date of purchase if the appliance is defective in workmanship or materials. To obtain service under this warranty, return the defective product to the service centre listed below, together with your purchase receipt. In the absence of a purchase receipt, the warranty period shall be 12 months from the date of manufacture.

ANY IMPLIED WARRANTIES, OBLIGATIONS OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED IN DURATION TO THE 12 MONTH DURATION OF THIS WARRANTY.

IN NO EVENT SHALL CONAIR BE LIABLE FOR ANY SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS OR IMPLIED, WHATSOEVER.

This warranty gives you specific legal rights, and you may also have other rights which vary from province to province.

CONAIR CONSUMER PRODUCTS ULC

SERVICE CENTRE:
Conair Consumer Products ULC
100 Conair Parkway,
Woodbridge, ON L4H 0L2

BM15C

For more information on any Conair Products call: 1-855-926-6247
Visit our Web site: www.conaircanada.ca
Consumer Call Centre E-mail: Consumer_Canada@Conair.com
©2021 Conair Consumer Products ULC
100 Conair Parkway, Woodbridge, ON L4H 0L2


IB-17314

21CN078956



BODY BENEFITS BY CONAIR®

ACUPRESSURE LUMBAR BELT WITH VIBRATION AND MASSAGE

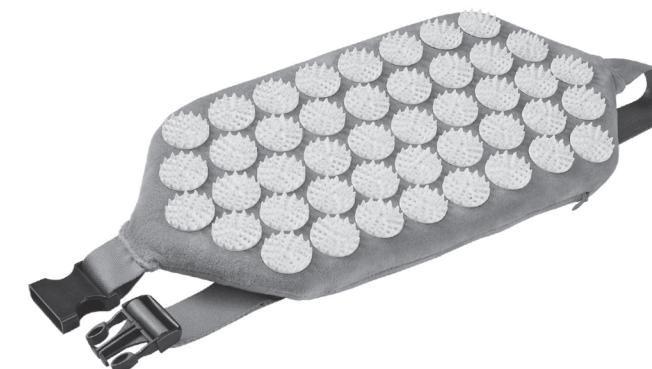

HAVE A QUESTION?

Please do not return this product to the retailer.

CALL US FIRST!

Our customers service and product experts are ready to answer **ALL** your questions.

Please call our **TOLL-FREE** customer service number at :

1-855-926-6247 Or visit us online at www.conaircanada.ca

Model BM15C Instruction Booklet

For your safety and continued enjoyment of this products,
always read the instruction book carefully before using.

IMPORTANT SAFETY INSTRUCTIONS

When using an electric massager, especially when children are present, basic precautions should always be taken, including the following:

READ ALL INSTRUCTIONS BEFORE USING THIS MASSAGER.

WARNING – To reduce the risk of burns, fire, electric shock, or injury to persons:

- 1) Do not use while bathing or in a shower.
- 2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 3) This massager is not intended for use by, on, or near children or people with certain disabilities.
- 4) Use this massager only for its intended use as described in this manual.
- 5) Keep the massager away from heated surfaces.
- 6) Never drop or insert any object into the massager.
- 7) To disconnect, turn all controls to the OFF position, then remove batteries from massager.
- 8) Do not use in bed.
- 9) Never use while sleeping.
- 10) **Never operate this massager if it is damaged, if it is not working properly, or if it has been dropped into water. Return the massager to an authorized service center for examination and repair.**
- 11) Do not use near any sharp objects to prevent massager from being punctured.
- 12) **CAUTION: CHOKING HAZARD.** This is not a toy. Not to be used as a baby or infant pillow. Do not open this massager for any reason. If this massager is damaged or opened, discard immediately, out of reach of children.
- 13) This massager contains a lumbar belt which has sharp plastic nodes on one side. Keep away from children and pets. Adults, use caution while handling.

SAVE THESE INSTRUCTIONS

About this massager

Your Lumbar Belt Massager helps the body to stimulate energy. It can be used to help ease tension and create deep relaxation. You may use the belt with or without the vibration feature.

OPERATING INSTRUCTIONS

Insert Batteries

Before using your massager you must insert 2 AA batteries (not included). Unzip cover to access the battery compartment. Insert batteries. Align positive (+) and negative (-) as indicated inside the battery compartment. Close zipper.

Place the massager on affected area to be treated. You can sit against the massager on a chair or sofa, but never sit directly on top of massager, as that could cause it to malfunction.

To activate the vibration massage, simply press the ON/OFF button on the back of the belt. (While the button is in the OFF position, you may use it without vibration.) You may use it on your lumbar back for up to 20 minutes. As a new user, you may want to use the lumbar belt over clothing. Once your body gets used to the nodes, you may use it on bare skin. Place the belt on your body with the nodes facing in, toward the body.

INSTRUCTIONS FOR USING BATTERIES

- Always purchase the correct size and grade of battery most suitable for the intended use.
- Replace all batteries of a set at the same time.
- Do not mix old and new batteries. Do not mix alkaline, standard (carbon-zinc), or rechargeable (nickel cadmium) batteries.
- Clean the battery contacts and also those of the device prior to battery installation.
- Ensure the batteries are installed correctly with regard to polarity (+ and -).
- Remove batteries from equipment that is not to be used for an extended period of time.
- Remove used batteries promptly.
- If you will not be using the massage feature for more than a few days, please remove batteries.

DISPOSAL OF BATTERIES

When replacing the batteries, make sure you dispose of the old batteries in accordance

with applicable regulations. They should be disposed of in a location specifically designated for that purpose so they can be recycled safely and without posing any risk for the environment. Do not burn or bury them.

CAUTIONS

- This massager should not be used in the following cases: Tuberculosis, benign and malignant tumors, hemorrhages, inflammations of the skin, phlebitis and thrombosis, open and fresh wounds, bruises, broken skin or varicose veins, undetermined calf pain or joint pain, or by persons with impaired sensitivity to heat.
- **Consult your physician should you have any questions regarding therapeutic use of this massager. Use of this product is solely at the user's risk and Conair Corporation assumes no liability for any loss or injury.**
- If you are pregnant, consult your physician before using.
- Do not use on any unexplained pain or swollen muscles, or following a serious muscular injury, before consulting your physician.
- Keep away from children.
- Discontinue use of massager if pain increases or swelling occurs after application.

MAINTENANCE INSTRUCTIONS

- 1) Store massager in a secure place away from excessive heat or sharp objects. If you will not be using the vibration feature for more than a few days, please remove batteries.
- 2) To clean the unit, remove the batteries and then wipe with a cloth. Do not use any strong chemical cleaners such as benzene.

RECOMMENDED MASSAGING TIMES

Generally, massage should not last more than 10 – 20 minutes on any one particular part of the body. Treatment may be applied once or twice daily as desired.

As individuals' metabolisms vary, so do their responses to massage. You will soon adjust the time to suit your particular needs and personal comfort. If you have any doubts regarding treatment, we recommend you consult your physician.