

## MOIST HEAT TREATMENT

1. Ensure heating pad is unplugged.
2. Wet the absorbent sponge and squeeze out excess water.
3. Insert the damp sponge between the heating pad and the cloth cover.
4. Plug the electrical cord into a 120V AC electrical outlet.
5. Push the ON/OFF button, the power indicator light will illuminate to signal that the heating pad is at Low setting. Press the settings buttons to choose among Low, Medium and High levels.
6. Place heating pad on top of the area to be treated.
7. When treatment is done, press the ON/OFF button to turn unit off and unplug. For maximum benefit, use should not exceed 30 minutes at one time.  
**Note:** This unit has a 60-minute auto-off feature.
8. After use, dry the cloth cover, absorbent sponge and the heating pad before storing.

**Note:** This heating pad has a 60-minute auto-off feature. During use, check skin under heating pad frequently. Do not leave heating pad unattended. After use, dry the cloth cover, absorbent sponge and the heating pad before storing.



## IMPORTANT SAFETY INSTRUCTIONS

### DO NOT DESTROY

**DANGER – To reduce the risk of burns, electric shock and fire, this product must be used in accordance with the following instructions:**

1. READ ALL INSTRUCTIONS.
2. **ALWAYS UNPLUG THE UNIT FROM THE ELECTRICAL OUTLET IMMEDIATELY AFTER USING.**
3. **DO NOT USE WHILE SLEEPING. AVOID FALLING ASLEEP WHILE USING THE PRODUCT.**
4. **DO NOT USE ON AN INFANT OR AN ANIMAL.**
5. **THIS PAD IS NOT TO BE USED ON OR BY ANY INDIVIDUAL WITH CERTAIN DISABILITIES, SLEEPING OR UNCONSCIOUS PERSON, OR A PERSON WITH POOR BLOOD CIRCULATION UNLESS CAREFULLY ATTENDED.**
6. **DO NOT USE ON AREAS OF INSENSITIVE SKIN.**
7. **BURNS CAN OCCUR REGARDLESS OF CONTROL SETTING; CHECK SKIN UNDER PAD FREQUENTLY.**
8. **NEVER USE PAD WITHOUT REMOVABLE COVER IN PLACE.**
9. Do not use in areas where oxygen is being administered.
10. **DO NOT USE PINS OR OTHER METALLIC MEANS TO FASTEN THIS PAD IN PLACE.**
11. **DO NOT SIT ON OR CRUSH PAD – AVOID SHARP FOLDS.**
12. **NEVER PULL THIS PAD BY THE SUPPLY CORD.**
13. **DO NOT USE THE CORD AS A HANDLE.**
14. CAREFULLY EXAMINE INNER COVER AND LINE CORD BEFORE EACH USE. DISCARD THE PAD IF INNER COVERING OR LINE CORD SHOW ANY

## RESTRICTED HEATING PAD USES

For additional warnings, please see the Important Safety Instructions and Caution sections.



Not for use by children



Do not sleep while using heating pad



Do not use heating pad bunched up



Do not use under covers



Do not sit on heating pad



Do not sharply fold or crease heating pad

### SIGN OF DETERIORATION.

15. Do not tamper with this pad in any way. There are no user serviceable parts. If, for any reason, this pad does not function satisfactorily, return to authorized service center at address given.
16. Do not use this pad with liniment, salve or ointment preparations that contain heat-producing ingredients, as skin burns could result. Use wherever hot applications are desirable for personal comfort and wherever recommended by your physician or dentist for the relief of pain.
17. **LOOP CORD LOOSELY WHEN STORING; TIGHT WRAPPING MAY DAMAGE CORD AND INTERNAL PARTS.**
18. Keep the pad away from heated surfaces.
19. Place pad on top and not under the part of the body needing heat.
20. Do not put this pad in the microwave.
21. Do not leave pad unattended, especially if children are present.
22. Do not use pads underneath bedcovers and other types of blankets. Do not place covers, such as sheets and blankets, above pads. Only the cloth cover provided with the pad should be used.
23. Do not use pad on sensitive skin or when taking medication that can cause sensitive skin.
24. Temperatures sufficiently high to cause burns may result from the improper use of this pad, regardless of the control switch setting; and therefore, if the pad is applied directly to the body, the patient's skin at the point of application should be inspected frequently.
25. Ensure that this pad is used only on a standard 120V circuit.
26. This heating pad must be kept dry and must not be operated in a wet or moist condition.
27. This pad is provided with a moisture-resistant envelope which should be carefully examined whenever the pad is put into service. If any signs of deterioration (such as checking, blistering, or cracking) are visible, the pad should not be used.
28. Never leave this pad connected to an electrical supply circuit when not in

actual use.

29. When not in use this pad should be kept in a cool, dark, dry place - preferably in the original carton.
30. This is not a toy. Do not allow children to touch and play it.

## SAVE THESE INSTRUCTIONS

### HEATING PAD COMPONENTS

- Machine Washable Cloth Cover (Figure 1)
- Heating Pad (Figure 2)
- Absorbent Sponge (Figure 3)



Figure 1

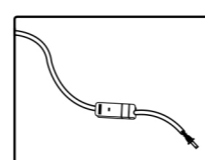


Figure 2

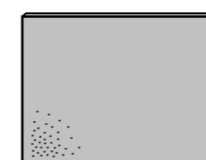


Figure 3

### Cautions:

- CAN ICES-3(B)/NMB-3(B)
- The cloth cover must always be used with the heating pad. An absorbent sponge is provided for moist heat application.
- This heating pad should not be used in the following cases: If you have diabetes, tuberculosis, benign or malignant tumors, phlebitis and thrombosis, hemorrhages, inflammation of the skin, undetermined calf pain or joint pain, or if you have poor circulation. Consult a physician before using if you have any kind of fungal infection, swelling, fracture or persistent pain.
- **Consult your physician should you have any therapeutic questions regarding use of this heating pad.**
- Do not use on any unexplained pain, or swollen muscles, or following a serious muscular injury before consulting your physician.
- **Keep away from children.**

## LIMITED THREE-YEAR WARRANTY (U.S. AND CANADA ONLY)

Conair will repair or replace, at our option, your appliance free of charge for 36 months from the date of purchase if the appliance is defective in workmanship or materials. To obtain service under this warranty, return the defective product to the service centre listed below, together with your sales slip and \$7.50 for postage and handling. In the absence of a receipt, the warranty period shall be 36 months from the date of manufacture. **ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, SHALL BE LIMITED IN DURATION TO THE 36-MONTH DURATION OF THIS WARRANTY.** **IN NO EVENT SHALL CONAIR BE LIABLE FOR ANY SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS OR IMPLIED, WHATSOEVER.** This warranty gives you specific legal rights, and you may also have other rights which vary from province to province.

### Conair Consumer Products ULC

**SERVICE CENTRE:**  
Conair Consumer Products ULC  
100 Conair Parkway,  
Woodbridge, ON  
L4H 0L2  
HP15RBWC

For more information on any Conair Products call:  
1-855-926-6247

Visit our Web site:  
[www.conaircanada.ca](http://www.conaircanada.ca)

Consumer Call Centre E-mail:  
[Consumer\\_Canada@Conair.com](mailto:Consumer_Canada@Conair.com)

©2018 Conair Consumer Products ULC

IB-15545

18CN033013

# BODY BENEFITS<sup>®</sup>

by **CONAIR<sup>®</sup>**



## MOIST / DRY KING SIZE HEATING PAD

### Have a Question?

*Please do not return this product to the retailer.*

### Call Us First!

**Our customer service and product experts are ready to answer ALL your questions.**

**Please call our TOLL-FREE customer service number at: 1-855-926-6247 or visit us online at [www.conaircanada.ca](http://www.conaircanada.ca)**

### Model HP15RBWC Instruction Booklet

For your safety and continued enjoyment of this products, always read the instruction book carefully before using.

## OPERATING INSTRUCTIONS

### DRY HEAT TREATMENT

This heating pad comes assembled (Figure 4).

This appliance has a polarized plug (one blade is wider than the other). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to defeat this safety feature. This heating pad is intended for household use only and it is designed to operate at 110-125 Volts AC. Use on Alternating Current (60 Hertz) only.

1. Plug the electrical cord into a 120V AC electrical outlet.
2. Push the ON/OFF button. The power indicator light will illuminate to signal that the heating pad is at Low setting. Press the settings buttons to choose among Low, Medium and High levels.
3. Choose your heat setting: LOW, MED, HIGH. Place heating pad on top of the area to be treated.
4. When treatment is done, press the ON/OFF button to turn unit off and unplug. For maximum benefit, use should not exceed 30 minutes at one time. **Note:** this unit has a 60-minute auto-off feature.
5. Do not use or store heating pad in a folded position. Use and store lying flat.

Burns can occur regardless of control setting; check skin under pad frequently.

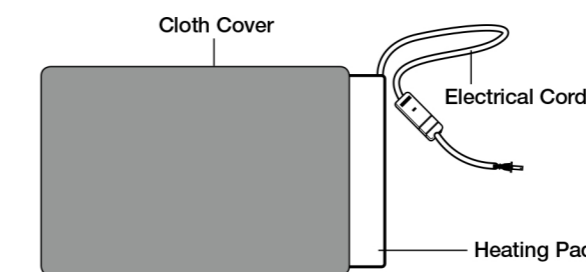


Figure 4

